


Single-Use Catheters for GI Manometry

Mediplus single-use gastrointestinal manometry catheters are ideal for all types of water-perfused manometry of the GI tract.

The range comprises single-channel, four-channel, eight-channel and 16-channel variants in external diameter sizes of 2.3mm, 3.2mm, 3.9mm, 4.5mm, 4.7mm and 4.9mm with a wide range of perfusion port configurations. All multi-lumen catheters have radiopaque markers at every eye (perfusion port) and black lines marked every centimetre to 85cm starting 1cm from the lumen nearest to the catheter tip.

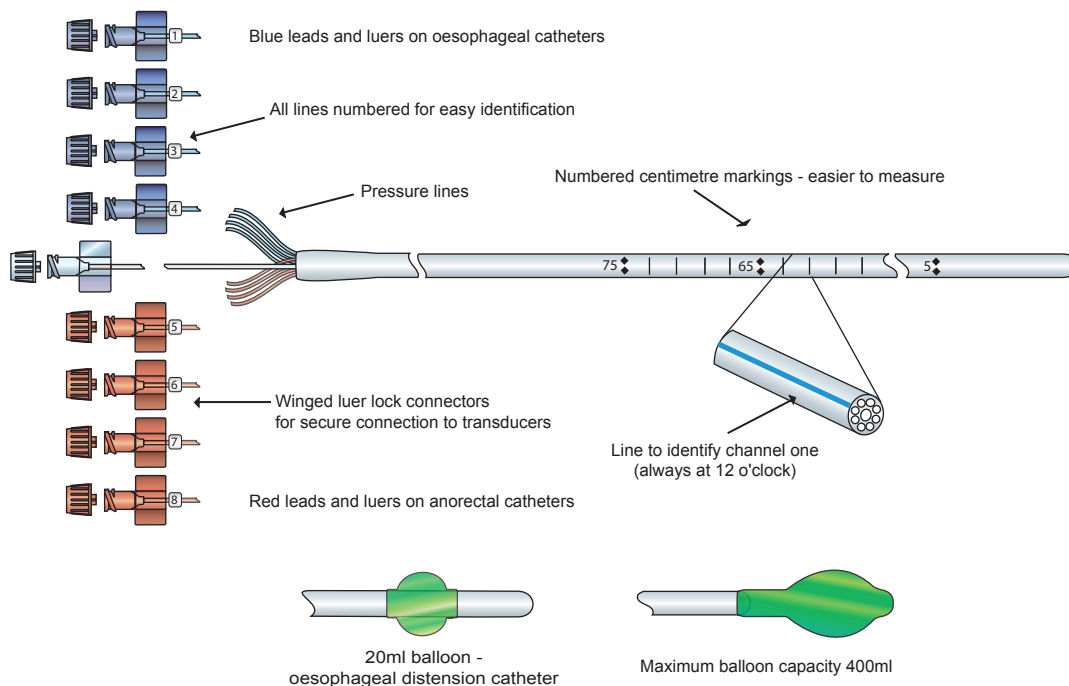

NOTE: Unless otherwise indicated, all catheters have open central lumen for guidewire insertion

Designed and priced for single-use - to remove risk of cross-infection, e.g. with *Helicobacter Pylori*, and exposure of staff to toxic chemicals used in cleaning reusable catheters, they offer the following benefits:

- ★ Small 2.3mm bore catheter - unique to Mediplus - ideal for paediatric applications
- ★ Multiple-channel measurement for a greater range of measurements possible than with microtips
- ★ Transparent tubing, so that any air bubbles in the line will be visible
- ★ Winged luer locks for secure connection to transducers
- ★ Numbered centimetre markings for easy measurement
- ★ Co-extruded stripe to locate channel one (except on anorectal catheter code 2304)

Patent No. GB2379389

Single-use Oesophageal Manometry Range

Code	Product Description	Configuration (cm)	Channels	Length	Tube o.d.	Pack Quantity
2244-1	Oesophageal Manometry	<p>5 5 5 5 5 cm Lumen # 8-5 4 3 2 1</p>	8	150cm	3.9mm	5
2244-1S	Oesophageal Manometry (Soft Catheter)	<p>5 5 5 5 5 cm Lumen # 8-5 4 3 2 1</p>	8	150cm	3.9mm	5
2244-8	Oesophageal Manometry	<p>5 5 5 5 5 cm Lumen # 1-4 5 6 7 8</p>	8	150cm	3.9mm	5
2280	Vector Manometry	<p>3 Lumen # 1-8</p>	8	150cm	3.9mm	5
2283	Oesophageal Distension Catheter (closed tip)	<p>5 20ml 5 5 5 5 cm Lumen # 8-5 4 3 2 1</p>	8	150cm	3.9mm	5
2243	Oesophageal Manometry	<p>5 5 5 5 5 cm Lumen # 8-5 4 3 2 1</p>	8	150cm	2.3mm	5
2240	Four Channel Oesophageal Manometry	<p>5 5 5 5 cm Lumen # 1 2 3 4</p>	4	180cm	2.3mm	5
2241	Four Channel Oesophageal Manometry	<p>3 3 3 3 cm Lumen # 1 2 3 4</p>	4	180cm	2.3mm	5

Patent No. GB2379389

Specialist Products in the Single-use Oesophageal Range

Mediplus Soft Oesophageal GI Catheter **Product code 2244-1S**

Mediplus Soft Oesophageal manometry catheters provide major benefits for both patients and user. The unique softer material formulation ensures easier intubation and **unlike** silicone the surface characteristics of our catheter material does not snag (catch) on gloves. The catheter is soft yet not flaccid, so there is no need for an introducer or stiffener. This could provide the Clinician with the facility of intubation without the need for an assistant in attendance.

Mediplus Paediatric Oesophageal GI Catheters **Product codes 2240, 2241 & 2243**

Mediplus Paediatric Oesophageal manometry catheters are ideal for all types of water perfused manometry of the GI tract. They offer a small 2.3mm outer diameter with eight or four channels.

Mediplus Oesophageal Distension Catheter **Product code 2283**

The small 3.9mm external diameter of the nine-lumen catheter allows for ease of intubation, whilst ensuring accurate pressure readings. A 20ml balloon can be inflated quickly and comfortably with the catheter in situ.

Please see previous page for product configurations and ordering details