


Double Lumen Cystometry

Mediplus double lumen cystometry catheters enable ease of patient catheterisation and simultaneous bladder filling and pressure measurement. The design is such that only 30% of the lumen cross-sectional area is devoted to cystometry measurement, with the remaining 70% dedicated to filling. This ensures accurate pressure measurements and a maximum fill rate of over 100 ml/minute*, with a catheter diameter smaller than would otherwise be needed; hence greater convenience and patient comfort.

The pressure lines on some catheters are 200 cm** in length to ensure direct connection of the patient's bladder to the pressure transducer without the need for additional extension tubing, so offering an accurate, cost-saving and convenient catheter. Suprapubic cystometry catheter set also available in Ch.6 and Ch.8. (see page U7 for details).

For greater user convenience we also offer Cystometry Sets comprising bladder pressure and filling catheters or bladder filling, pressure and rectal catheters in one pack. See page U6 for details.

*Ch.8 and Ch.10 only
** Nominal

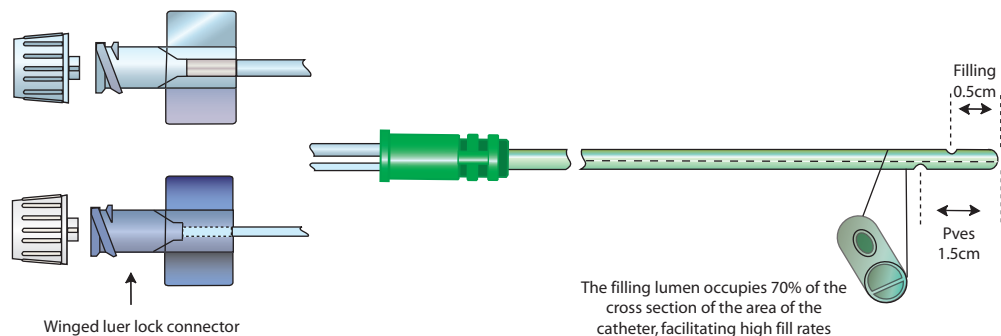

Note: Product Codes 5732, 5736 & 5738 have marker lines at 1cm intervals to 25cm

Code	Product Description	Lengths			Pack Quantity
Double Lumen Cystometry (Medical grade PVC - suitable for patients with latex hypersensitivity)					
		Filling Line	Pressure Line	Multi Lumen Length	
5716	6 Fr./Ch.	15cm	130cm	65cm	10
5717	6 Fr./Ch. (short)	15cm	130cm	20cm	10
5702	8 Fr./Ch.	15cm	130cm	65cm	10
5703	8 Fr./Ch. (short)	15cm	130cm	30cm	10
5712	10 Fr./Ch.	15cm	130cm	65cm	10
5713	10 Fr./Ch. (short)	15cm	15cm	65cm	10
5732	8 Fr./Ch. (standard)	7cm	7cm	30cm	10
5736	6 Fr./Ch. (standard)	7cm	7cm	30cm	10
5738	10 Fr./Ch. (standard)	7cm	7cm	30cm	10

These products are supplied sterile.