


Suprapubic Foley Catheter Introducing Set

Using the Seldinger catheter introduction system to place the Foley urinary drainage catheter offers a number of advantages to both the clinician and the patient:

- Less pressure needed to insert the trocar, so more comfortable under local anaesthetic and less stressful for the surgeon; high degree of control for accurate placement.
- Minimum trauma for the patient and minimum tissue damage.
- Locates the bladder with 18g hypodermic needle, so more confidence of inserting the trocar into the bladder, as the track has already been secured by the guide wire. - No more guessing the angle the anaesthetic needle is pointing in!
- Three-stage guide wire improves insertion and removal and guarantees insertion of the trocar along the anaesthetised track. The variable stiffness of the guide wire allows secure and comfortable curling up in the bladder while reliably guiding the trocar and preventing damage to the posterior wall.
- The soft, uniformly inflated balloon is flush to the catheter shaft to improve insertion and removal, and sits well against the bladder.

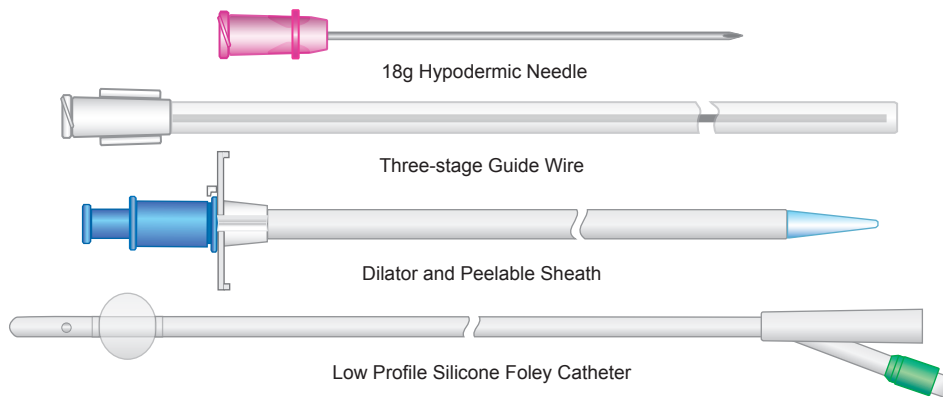

Set includes 2 x 10ml syringe and a scalpel

Code	Product Description	Fr./Ch.	Length	Pack Quantity
5754	Silicone Suprapubic Foley Catheter Introducing Set (<i>Latex Free</i>)	14	43cm	5
5756	Silicone Suprapubic Foley Catheter Introducing Set (<i>Latex Free</i>)	16	43cm	5

Replacement Catheters

5764	Silicone 2-way Foley Catheter	14	43cm	10
5766	Silicone 2-way Foley Catheter	16	43cm	10
5768	Silicone 2-way Foley Catheter	18	43cm	10