


UPP: Triple Lumen Catheters

This three-lumen design enables cystometry and UPP with simultaneous perfusion to be performed. The three-lumen section of the catheter has a urethral pressure profile cross sectional area showing 50% apportioned to filling for improved filling rates, 25% to cystometry and 25% to UPP. The two colour-coded Ch.4.5 pressure lines (red and blue on the Ch.9 catheters and yellow and blue on the Ch.7 catheter) are an additional 130cm long to enable direct bladder-transducer connection without the need for extension tubing.

The PVC Ch.9 triple lumen catheter is available in three variants, codes 5331 and 5334 with 130cm pressure lines, and code 5332 with shorter 12cm pressure lines.

The Polyurethane Ch.7 triple lumen catheter, code 5333, has two collars by the UPP port to enable a reading of the average pressure through 360° to be obtained.

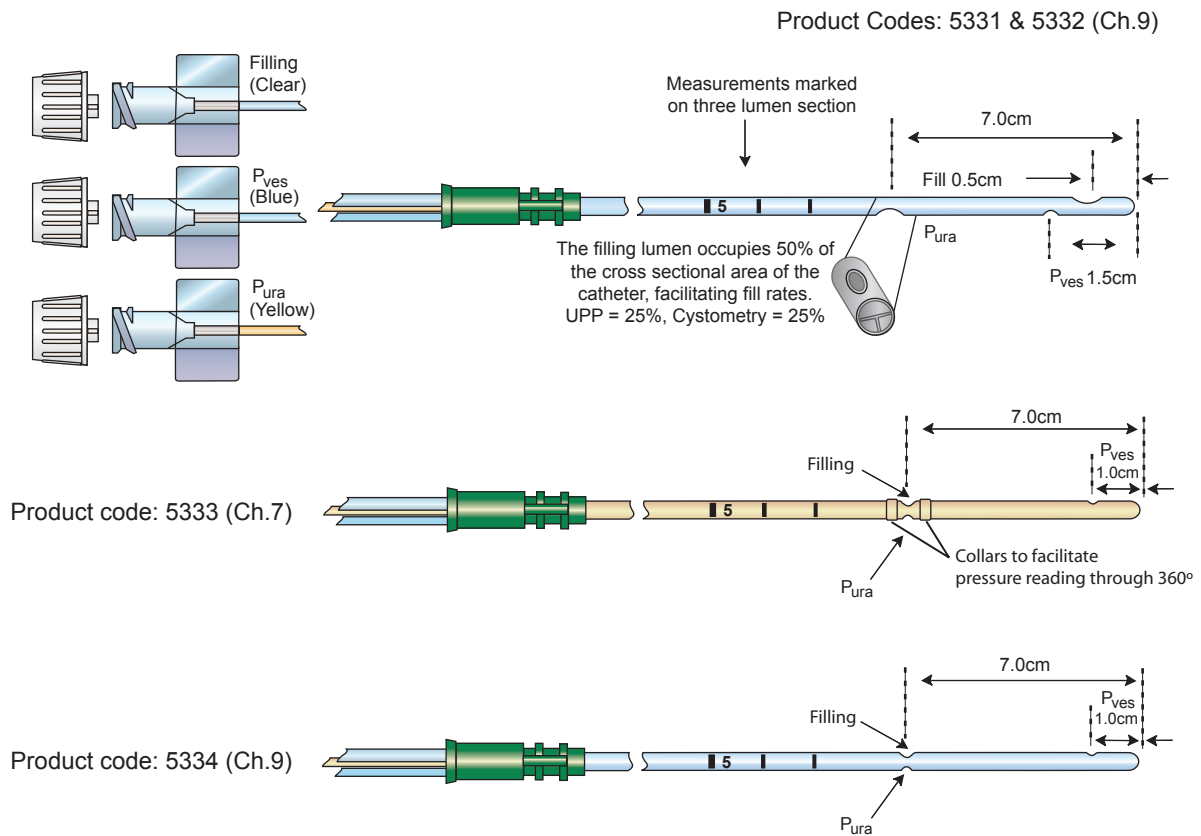

Code	Product Description	Lengths	Pack Quantity
Triple Lumen UPP Catheters - Sterile (Medical grade PVC and Polyurethane - suitable for patients with latex hypersensitivity)			
		Filling Line	P_{ves} & P_{ura} Line
		Multi Lumen Length	
5331	9 Fr./Ch. PVC	15cm	130cm
5332	9 Fr./Ch. PVC	15cm	12cm
5333	7 Fr./Ch. Polyurethane	12cm	130cm
5334	9 Fr./Ch. PVC	15cm	130cm

These products are supplied sterile.